



VELA[®] ventilator


Comprehensive. Versatile. Proven.

Healthcare providers need solutions that empower them to deliver quality patient care. The VELA ventilator empowers clinicians with a comprehensive range of modalities for pediatric and adult patients.



Advancing ventilation technology with an all-in-one solution

In busy hospitals where versatility is critical, the VELA ventilator empowers clinicians with a comprehensive range of modalities for pediatric and adult patients requiring either invasive or noninvasive ventilatory support. Whether in the Emergency Room, ICU or Long Term Care the VELA ventilator has the performance capability and versatility to manage a broad range of patient care settings throughout the hospital.

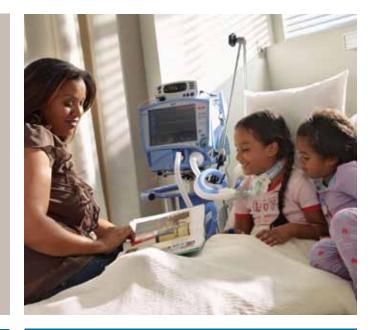

Comprehensive

- Invasive and noninvasive critical care ventilator
- Pediatric through adult applications
- Advanced settings to treat the most difficult to ventilate patient
- Multiple noninvasive modes with automatic leak compensation to provide more ventilatory support options
- Dual limb noninvasive ventilation helps eliminate CO₂ and rebreathing of exhaled gases
- Apnea back up ventilation with user selectable breath type
- Intuitive user interface facilitates ease of use
- Integrated end tidal capnography monitoring

Versatile

- Performs multiple applications from the Emergency Room to Critical Care and sub-acute settings
- Operates independent of air sources for greater flexibility
- Portable unit allows patients to stay on NIV therapy during transport
- High pressure oxygen inlet with blender
- Low pressure oxygen inlet with accumulator
- Internal battery capacity of six hours
- Compatible with speaking valves

Proven

- Patented turbine technology with flow capabilities of 180 LPM
- Two-year warranty, with a five-year or 40,000-hour warranty on the turbine
- Accurate, integrated FiO₂ delivery, monitoring and alarms
- Complete 24 hour trended data of all patient parameters
- No scheduled maintenance of the turbine with a low cost of ownership

Intuitive operation

The Simple Touch[™] user interface in the VELA ventilator reduces training time, produces consistency in method and enhances clinician confidence.

Mode ~

Visual display of current mode with easy and flexible access to all modes of operation

User selectable patient s monitor values

Five monitored parameters selectable from 19 options are displayed continuously on main screen

Main screen

An active matrix, full color LCD touch screen features scaleable graphic displays

Ventilation settings


Ventilator settings are displayed on all screens

Nebulizer

Synchronized on inspiration, nebulization time can be set in increments of one minute for up to a maximum of 60 minutes

Main Access additional menus easily from the main screen Alarm message window Alarm messages are displayed

visually with audible sound when active according to priority

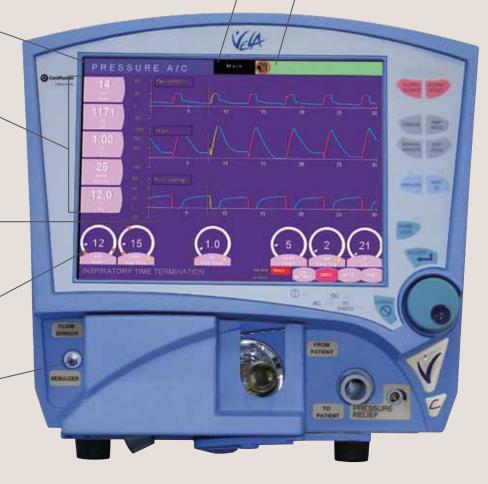



Complete trended data

An intelligent, scalable display organizes minute-byminute data into trended histograms and tables for all monitored parameters.



Dedicated monitoring screen

Configurable monitoring screen allows you to match your critical care documentation.



Advanced monitoring options

Integrated end tidal capnography provides clinicians with advanced monitoring to evaluate therapy and assess respiratory disease and airway placement. \triangle WARNING: U.S. Federal Law restricts this device to sale by or on the order of a physician.

CareFusion 22745 Savi Ranch Parkway Yorba Linda, CA 92887

> 800.231.2466 toll-free 714.283.2228 tel 714.283.8493 fax

C € ••••

CareFusion Germany 234 GmbH Leibnizstrasse 7 97204 Hoechberg Germany +49 931 4972-0 tel

+49 931 4972-423 fax

© 2010 CareFusion Corporation or one of its subsidiaries. All rights reserved. VELA and Simple Touch are trademarks or registered trademarks of CareFusion Corporation or one of its subsidiaries. RC1021 (0910/4000) L 3145 Rev. B

CareFusion Yorba Linda, CA

carefusion.com

